


Microsphere Selection

Suggestions and strategies for choosing the best beads for common microsphere applications.

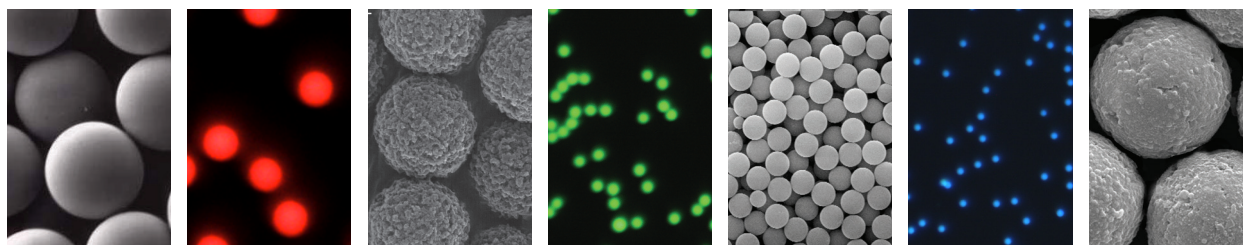

Microspheres offer a highly convenient and flexible system for developing reagents for assays and bioseparations, and for use as instrument standards. With many varieties of microspheres available, it is important to think about the demands the application will place on them when choosing a base bead. Physical and optical properties should be considered in the context of handling and detection, and thought should also be given to requirements for diameter and size distribution, composition, surface chemistry, and any other needed properties.

Diameter

Microsphere size may be critical to the proper function of an assay, or it may be a secondary consideration. Considering traditional diagnostic methods, the test or assay format commonly dictates particle size, such as the use of submicron spheres ($\leq 0.5\mu\text{m}$) in lateral flow tests and turbidimetric assays, or the use of cell-sized spheres ($\sim 4\text{-}10\mu\text{m}$) for bead-based flow cytometric assays. In magnetic separations, the exact size of the magnetic particle may be unimportant provided that the particles are in some general size range, and offer desired separation characteristics.

Diameter also determines surface area. Small-diameter spheres present more surface area per unit weight, while larger spheres present more surface area per bead. Size also impacts ease of handling, processing considerations, and the amount of reagent needed for coating.

Composition

Common microsphere compositions include polystyrene (PS), poly(methyl methacrylate) (PMMA), and silica. These materials possess different physical and optical properties, which may present advantages or limitations for different applications.

Polymer beads are generally hydrophobic, and as such, have high protein binding abilities. However, they often require the use of some surfactant in the storage buffer to ensure ease of handling. During synthesis, functional monomers may be co-polymerized with styrene or methyl methacrylate to develop beads with surface reactive groups. Functional groups may be used in covalent binding reactions, and aid in stabilizing the suspension.

Silica microspheres are inherently hydrophilic and negatively charged. Consequently, aqueous silica suspensions rarely require use of surfactants or other stabilizers. Carboxyl- and amine-functionalized silica microspheres are available for use in common covalent coating protocols, and plain silica spheres may be modified using a variety of silanes to generate functional groups or alter surface properties.

Coating

Microspheres may be coated with ligands, such as antibodies or oligonucleotides, for use in diagnostic or separation applications. Microsphere coatings are typically optimized to achieve desired specific activity, while minimizing nonspecific interactions. Consideration should be given to the required stability, development time frame and budget, and the specific biomolecule to be coated. These factors will help determine the most fitting coating strategy (adsorption, covalent coupling, or affinity binding) for short- and long-term objectives.

Adsorption

Adsorption relies primarily on hydrophobic interactions between the biomolecule and polymer particle. Such coatings are fairly simple to conduct, involving incubation of the microspheres with the purified biomolecule. They typically require little optimization, and reagents may be developed relatively quickly. However, as adsorption relies on the formation of multiple attachment points between the molecule and particle, this strategy is typically reserved for use with proteins and non-functionalized polymer spheres. Adsorption is generally not suitable for hormones, peptides, or nucleic acids in hybridization-based applications, and protein adsorption to silica is expected to be less efficient than to polymer.

Covalent Coupling

Covalent coupling results in the permanent attachment of the molecule to the functionalized microsphere. It can provide needed stability when developing a commercial reagent, and for multiplexed assays, where analyte-specific bead populations are mixed. Additionally, specialized chemical linkers may be employed to address steric effects or to optimally orient the molecule. Although covalent binding protocols often involve a higher level of optimization than other approaches, coupling kits are available to simplify the process.

Affinity Binding

Affinity binding is a straightforward method for immobilizing primary antibodies or biotinylated molecules. Proteins A and G and Fc-specific antibody coatings permit the directed immobilization of primary antibodies, and streptavidin is used extensively for the binding of biotinylated molecules, such as antibodies, peptides, and oligonucleotides.

It is important to note that each binding strategy has benefits and limitations, which should be weighed in the context of study objectives and the demands that will be placed on the finished reagent.

Special Properties

Many applications in the life sciences demand added properties, such as fluorescence or a visible color, or iron oxide inclusions for magnetic separations. Polymer spheres (and polymer-based magnetic spheres) are often internally dyed via organic solvent swelling, and many standard products are available. Dye concentrations can be adjusted to produce beads with different intensities to meet special needs, such as QuantumPlex™ for multiplexed flow cytometric assays, or our Dragon Green or Flash Red Intensity Standards, which support imaging applications and associated instrument QC. Many surface- or internally-labeled fluorescent beads are also available as specialized flow cytometry standards.

Various types of superparamagnetic microparticles are available as well – with different matrices, magnetite content, surface groups, etc. For new assays or applications, magnetic beads should be evaluated with application demands in mind.

For additional information regarding microsphere selection, please see our **Microsphere Reagent Development Booklet** or visit our website at www.bangslabs.com.





Bangs Laboratories manufactures a variety of magnetic, polymeric and silica microsphere products setting the standards for diagnostic, research, and flow cytometry applications. No matter the project, we have a product that serves or we'll work to custom-design a solution to fit. And that's not the half of it.


We also stand behind our products. Regardless of the size of your question or the size of your company, we offer tech support, absolutely free.

Sound interesting? 

Visit: www.bangslabs.com

 @particledoc

 info@bangslabs.com

 800.387.0672